


AcceX

0.014" PTCA Guidewire


The AcceX PTCA Guide wire Series are intended to facilitate the placement of devices during diagnostic and interventional procedures. The guide wires have a distal section made of hydrophilic coated coil, the proximal section is made with a SS hardened core wire. The core wire is fixed at the proximal end and extends to tip of the distal end. The proximal section of the corewire is PTFE coated. The distal segment of the core wire has a profiled taper. The core wire enhances the stiffness of the guide wire. A coil is provided at the distal section of the core wire up-to 30 cm, which shall provide enough flexibility to the distal section of the guidewire. The tip of the guidewire is constructed with a Platinum alloy marker for visualization under fluoroscopy.

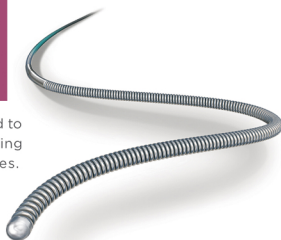

The AcceX PTCA guide wire series are offered in 0.014 inch outer diameter with different tip configuration (straight or J tip) with various support of Floppy, Regular and extra support and are available in length 185cm.

INTENDED USE OF DEVICE

The PTCA Guide Wire is used to facilitate the placement of devices during diagnostic and interventional procedures. The guide wires are intended to facilitate catheter placement and exchange during diagnostic or interventional procedures, where increased support, distal flexibility and low surface friction of the guide wire is needed.

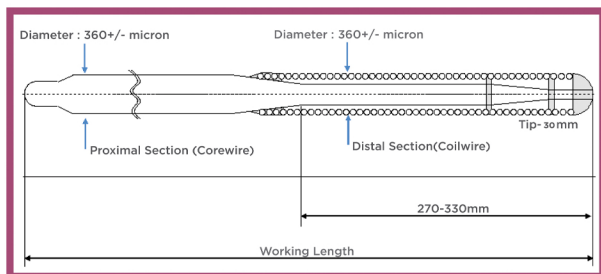

KEY PARAMETERS

- Unicore design
- Multiple core support options
- Hydrophilic Distal shaft
- Excellent Maneuverability
- Good Support for device delivery
- Excellent trackability
- Good Workhorse wire

Proximal Shaft	Stainless Steel
Distal Shaft	Stainless Steel
Tip	Platinum Alloy
Coating - Proximal	PTFE coating
Coating - Distal	Hydrophilic coating

Guide Wire Variants	Wire Diameter	Working Length	Tip Shape	Tip Load	Tip Length	Vasmed Code
Accex Floppy	0.014"	185 cm	J tip/Straight tip	<0.6g	3 cm	VM-AX-014-FJ-185 VM-AX-014-FS-185
Accex Regular	0.014"	185 cm	J tip/Straight tip	0.7g to 0.9g	3 cm	VM-AX-014-RJ-185 VM-AX-014-RS-185
Accex Extra Support	0.014"	185 cm	J tip/Straight tip	>0.9g	3 cm	VM-AX-014-ES-185 VM-AX-014-EES-185